## MultiMix™

### Triple-Colour Reagent

#### Anti-Human Kappa Light Chains/FITC

#### Anti-Human Lambda Light Chains/RPE

#### Anti-Human CD19/RPE-Cy5

#### Code No./ Code/ Code-Nr. TC051

ENGLISH
---------

#### Intended use

For in vitro diagnostic use.

TC051 is intended for use in flow cytometry for simultaneous detection and enumeration of kappa light chains and lambda light chains on human B cells. In flow cytometry, antibodies to kappa light chains are useful for the demonstration of cell surface kappa light chains, and antibodies to lambda light chains are useful for the demonstration of cell surface lambda light chains, and, thus, for the identification of monoclonality (clonal excess) in B-cell lymphoproliferative disorders together with a panel of other antibodies (1-4). CD19 is the broadest lineage-specific surface marker for B cells and it is present on the surface of virtually all B lymphocytes, including early B progenitor cells (5). CD19 expression is maintained in B-lineage cells that have undergone neoplastic transformation (6). Antibodies to CD19 are considered essential for the initial evaluation of acute and chronic lymphoproliferative disorders (7). Interpretation of results must be made within the context of the patient's clinical history and other diagnostic tests by a certified professional.

#### Reagent provided

TC051 comprises the following three, carefully matched, fluorescent antibodies:

Polyclonal Rabbit Anti-Human Kappa Light Chains conjugated with fluorescein isothiocyanate isomer 1 (FITC).

Polyclonal Rabbit Anti-Human Lambda Light Chains, conjugated with R-phycoerythrin (RPE).

Monoclonal Mouse Anti-Human CD19, Clone HD37, conjugated with conjugated R-phycoerythrin-Cyanine 5 (RPE-Cy5).

The Anti-Kappa Light Chains conjugate has been produced from a F(ab')<sub>2</sub> fragment of affinity-isolated polyclonal rabbit antibody. The Anti-Lambda Light Chains conjugate has been produced from a F(ab')<sub>2</sub> fragment of affinity-isolated polyclonal rabbit antibody. The Anti-CD19 conjugate has been produced from a purified monoclonal mouse antibody. TC051 is provided in liquid form in buffer containing 1% bovine serum albumin (BSA) and 15 mmol/L NaN<sub>3</sub>, pH 7.2. Each vial contains 50 tests (10 µL of conjugate for up to 10<sup>6</sup> leucocytes from normal human peripheral blood).

Antibody Code No.	Fluorochrome	Control Reagent Code No.
TC051	FITC, RPE and RPE-Cy5	X0957

#### Immunogens

Rabbit Anti-Human Kappa Light Chains: Polyclonal immunoglobulin light chains of kappa type isolated from a pool of human sera.

Rabbit Anti-Human Lambda Light Chains: Polyclonal immunoglobulin light chains of lambda type isolated from a pool of human sera.

Anti-CD19, Clone HD37: Hairy cell leukaemia cells.

#### Specificity

**Anti-Kappa Light Chains:** Reacts with free kappa chains as well as kappa chains in intact immunoglobulin molecules. The specificity has been ascertained as described below. To obtain maximum sensitivity, the crossed immunoelectrophoresis specificity test was performed prior to affinity purification and pepsin degradation.

**Crossed immunoelectrophoresis:** Only kappa-related precipitates appear when the antibody is run against human plasma. Staining: Coomassie Brilliant Blue.

**Flow cytometry:** When Anti-Kappa Light Chains/FITC is applied as described in the staining procedure in combination with Anti-CD19/RPE, HD37, on lysed human whole blood, a specific staining of a part of the CD19-positive B lymphocytes is seen corresponding to the expected range of kappa light chain expression.

**Anti-Lambda Light Chains:** Reacts with free lambda chains as well as lambda chains in intact immunoglobulin molecules. The specificity has been ascertained as described below. To obtain maximum sensitivity, the crossed immunoelectrophoresis specificity test was performed prior to affinity purification and pepsin degradation.

**Crossed immunoelectrophoresis:** Only lambda-related precipitates appear when the antibody is run against human plasma. Staining: Coomassie Brilliant Blue.

**Flow cytometry:** When Anti-Lambda Light Chains/RPE is applied as described in the staining procedure in combination with Anti-CD19/FITC, HD37, on lysed human whole blood, a specific staining of a part of the CD19-positive B lymphocytes is seen corresponding to the expected range of lambda light chain expression.

Flow cytometric analysis of single suspensions from formalin-fixed, paraffin-embedded tissue samples demonstrated that Anti-Lambda Light Chains labels reactive hyperplastic lymph nodes, and so does Anti-Kappa Light Chains (10/10 cases). In non-Hodgkin's lymphomas, the antibody labeled 5/10 cases. The remaining cases were positive for kappa light chains (4/10 cases) or showed no expression of light chains (1/10 cases) (4).

Flow cytometric analysis of peripheral blood lymphocytes demonstrated that Anti-Kappa Light Chains and Anti-Lambda Light Chains, respectively, label a proportion of B-cell chronic lymphocytic leukaemias. Thus, in one study of 121 cases (2), 37 were positive for kappa and 17 were positive for lambda. In another study of 165 cases (3), 97 were positive for kappa and 68 were positive for lambda.

**Anti-CD19, HD37:** Was included in the Second, Third, Fourth and Fifth International Workshops and Conferences on Human Leucocyte Differentiation Antigens, and studies by a number of laboratories confirmed its reactivity with CD19 (8). Anti-CD19, HD37, labels human B cells in peripheral blood, bone marrow and other tissues. Additionally, B-cell lymphoproliferative disorders gave positive reactions with the antibody, i.e. acute lymphoblastic leukaemia, chronic lymphocytic leukaemia, hairy cell leukaemia, lymphoblastic lymphoma (Burkitt type), centroblastic/centrocytic lymphoma, and diffuse non-Hodgkin's lymphoma. The antibody was shown to be unreactive with other cells in the human haematopoietic system and did not react with non-haematopoietic cells, e.g. in kidney, liver, breast or lung tissues (8).

**TC051, Flow cytometry:** When TC051 is applied as described in the staining procedure, on lysed human whole blood, a specific staining of a part of the CD19 positive B lymphocytes is seen, within the expected range of kappa and lambda light chains expression.

#### Precautions

- For professional users.
- This product contains sodium azide (NaN<sub>3</sub>), a chemical highly toxic in pure form. At product concentrations, though not classified as hazardous, sodium azide may react with lead and copper plumbing to form highly explosive build-ups of metal azides. Upon disposal, flush with large volumes of water to prevent metal azide build-up in plumbing.
- As with any product derived from biological sources, proper handling procedures should be used.

#### Storage

Store in the dark at 2-8 °C. Do not use after expiration date stamped on vial. If reagents are stored under any conditions other than those specified, the conditions must be verified by the user. There are no obvious signs to indicate instability of this product. Therefore, positive and negative controls should be run simultaneously with patient specimens. If unexpected staining is observed which cannot be explained by variations in laboratory procedures and a problem with the antibody is suspected, contact our Technical Services.

#### Product-specific limitations

It has been observed that RPE-Cy5-conjugates may bind to monocytes, resulting in background staining (9).

#### Staining procedure

Before staining samples of peripheral blood, the mononuclear cells must be isolated by centrifugation on a separation medium or the blood sample washed to remove the soluble serum proteins. As human monocytes bind serum immunoglobulins via their surface Fc receptors, these cells should be removed by depletion or identified.

- Transfer 100 µL of anticoagulated (EDTA) blood to a 12 mm x 75 mm polystyrene test tube.
- Add 2 mL of PBS (DakoCytomation code No. S3024). Mix gently by using a vortex mixer.
- Centrifuge the tube at 300 x g for 5 minutes. Gently aspirate the supernatant and discard it leaving approximately 50 µL of fluid in the tube.
- Repeat steps 2 and 3 two more times.
- Add 10 µL Rabbit Immunoglobulin Fraction (DakoCytomation code No. X0903) for blocking. Mix gently by using a vortex mixer and incubate the tube at 37 °C for 30 minutes.
- Add 10 µL of TC051 and mix gently by using a vortex mixer.
- Incubate the tube in the dark at 2-8 °C for 30 minutes or at room temperature (20-25 °C) for 15-30 minutes.
- Add 100 µL of Uti-Lyse™ Reagent A (DakoCytomation code No. S3325) to the tube and mix gently by using a vortex mixer. Incubate the tube for 10 minutes at room temperature in the dark.
- Add 1 mL of Uti-Lyse™ Reagent B (DakoCytomation code No. S3325) to the tube and mix gently by using a vortex mixer. Incubate the tube for 10 minutes at room temperature in the dark.
- Centrifuge the tube at 300 x g for 5 minutes. Gently aspirate the supernatant and discard it leaving approximately 50 µL of fluid in the tube.
- Add 2 mL of PBS to the tube and resuspend the cells by using a vortex mixer.
- Repeat step 10.
- Resuspend the cells in an appropriate fluid for flow cytometry, e.g. 0.3 mL PBS.
- Analyse the sample on a flow cytometer or store it in the dark at 2-8 °C until analysis. Samples should be analysed within 24 hours after staining.

Note that fluorochrome conjugates are light sensitive, and samples should be protected from light during the staining procedure and until the analysis.

#### Procedural notes

Step 1: It is optional to include a suitable positive and negative control sample with each run for reagent and preparation control.

Step 6: The volume of conjugate recommended is a guideline only. Optimal staining conditions may vary depending on specimen and preparation method and should be determined by each individual laboratory.

It is optional to include a control reagent test tube. The control reagent should match the conjugated antibody reagent. The recommended Triple-Colour Control Reagent for TC051 is DakoCytomation code No. X0957.

Steps 8 and 9: If another cell-lysing reagent is used, please follow the recommendations for that reagent. Note that if the alternative lysing reagent does not contain fixative, such as DakoCytomation EasyLyse™, code No. S2364, the PBS in step 13 should contain 1% paraformaldehyde unless the sample is analysed within the time frame recommended for the lysing reagent.

Step 14: Abnormal numbers of cells expressing the target antigen or cells showing aberrant expression levels of the antigen can be expected in some disease states. This may result in an altered pattern of staining. It is important to understand the normal expression pattern of an antigen and its relationship to the expression of other relevant antigens in order to perform appropriate analysis.

Multicolour reagents are preferable to single-colour reagents for the comprehensive analysis of flow cytometry specimens. The correct use of colour compensation is particularly important in multicolour analysis.

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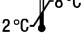
1. Johnson A, Olofsson T. Flow cytometric clonal excess analysis of peripheral blood, routine handling, and pitfalls in interpretation. *Cytometry* 1993;14:188-95.
2. Cartron G, Linassier C, Bremond JL, Desablens B, Georget MT, Fimbel B, et al. CD5 negative B-cell chronic lymphocytic leukemia: clinical and biological features of 42 cases. *Leuk Lymphoma* 1998;31:209-16.
3. Gandini D, Lanza F, Latorraca A, Levato F, Del Senno L, Castoldi G. Immunophenotypic and genotypic characterization of B-cell chronic lymphocytic leukemia patients from Northern Italy. *Haematologica* 1993;78:18-24.
4. Leers MPG, Theunissen PHMH, Ramaekers FCS, Schutte B, Nap M. Clonality assessment of lymphoproliferative disorders by multiparameter flow cytometry of paraffin-embedded tissue: an additional diagnostic tool in surgical pathology. *Hum Pathol* 2000;31:422-7.
5. Pezzutto A, Dörken B, Rabinovitch PS, Ledbetter JA, Moldenhauer G, Clark EA. CD19 monoclonal antibody HD37 inhibits anti-immunoglobulin-induced B cell activation and proliferation. *J Immunol* 1987;138:2793-9.
6. Sato S, Tedder TF. BC3. CD19 workshop panel report. In: Kishimoto T, Kikutani H, von dem Borne AEG, Goyert SM, Mason DY, Miyasaka M, et al., editors. *Leucocyte typing VI. White cell differentiation antigens. Proceedings of the 6th International Workshop and Conference*; 1996 Nov 10-14; Kobe, Japan. New York, London: Garland Publishing Inc.; 1997. p. 133-5.
7. Braylan RC, Orfao A, Borowitz MJ, Davis BH. Optimal number of reagents required to evaluate hematolymphoid neoplasias: results of an international consensus meeting. *Cytometry* 2001;46:23-7.
8. Mason DY, Ladyman H, Gatter KC. Immunohistochemical analysis of monoclonal anti-B cell antibodies. In: Reinherz EL, Haynes BF, Nadler LM, Bernstein ID, editors. *Leukocyte typing II. Proceedings of the 2nd International Workshop on Human Leukocyte Differentiation Antigens*; 1984 Sept 17-20; Boston, USA. New York, Berlin, Heidelberg, Tokyo: Springer-Verlag; 1986. Volume 2. p. 245-55.
9. van Vugt MJ, van den Herik-Oudijk IE, van de Winkle JGJ. Binding of PE-CY5 conjugates to the human high-affinity receptor for IgG (CD64). *Blood* 1996;88:2358-61.

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